



Five Essential Tools
to Achieve Surgical Convenience

Assist Kit

- 1. Tissue Punch
- 2. Bone Profiler
- 3. Trepine Drill
- 4. Crest Remover
- 5. Autobone Collector




HIOSSEN
IMPLANT

1 Tissue Punch



Minimizes incisions and reduces healing time

- The guide pin configuration enables precise tissue removal at the designated location
- Recommend using a diameter 0.7-1.5mm smaller than the healing abutment diameter

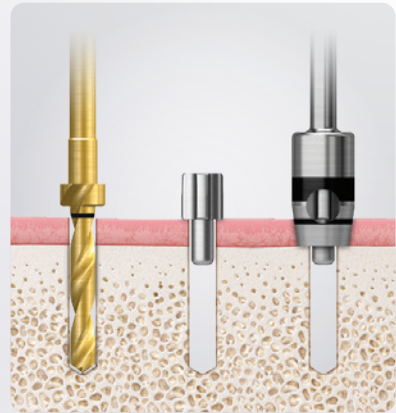


2 Bone Profiler




Removes protruding bone without internal damage of the implant

- Easily removes adjacent obstructive bone around the implant
- Detachable guide protects interior of the implant while removing the bones


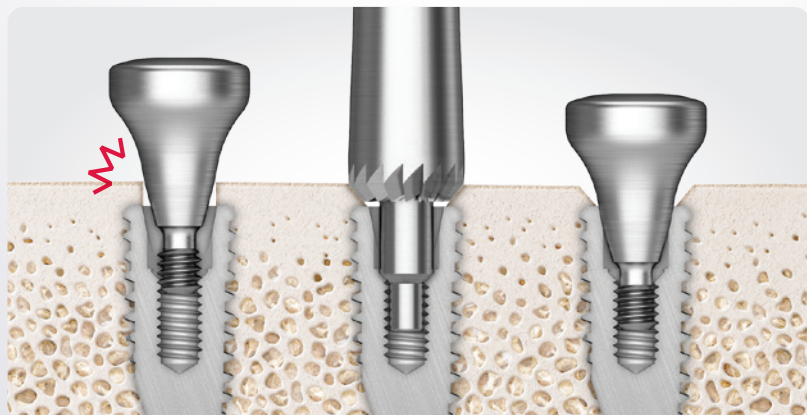


3 Trepphine Drill

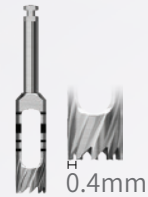

Thin blade wall (0.4mm) minimizes amount of the bone removed

- The drill is designed to harvest larger block bone from the donor site or to remove the failed implant
- The 0.4mm thin wall design reduces the quantity of removed bone

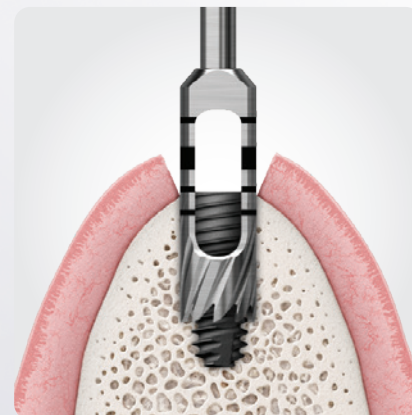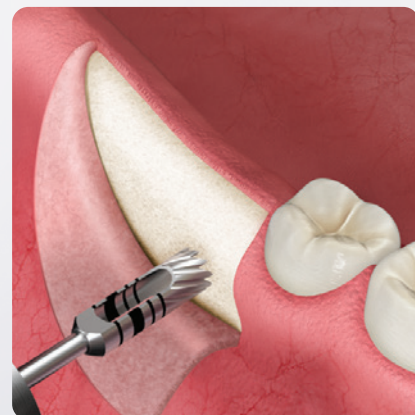


4 Crest Remover



Net-shaped blades of the bur speeds up ridge widening

- The blade is designed to minimize bone splattering effectively during the drilling process

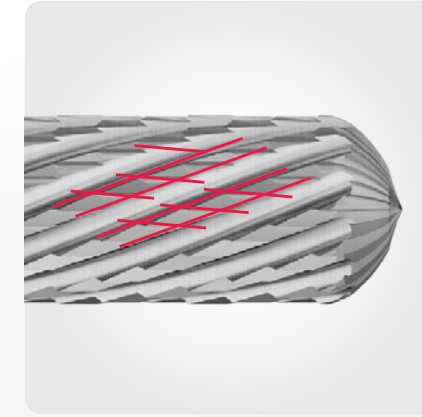


5 Autobone Collector



Collects bone chips through drilling

- Efficiently gather autogenous bone chips with a strong cutting force and design
- The stopper enables a safer method for collecting bone material (0.18g acquired at one operation)

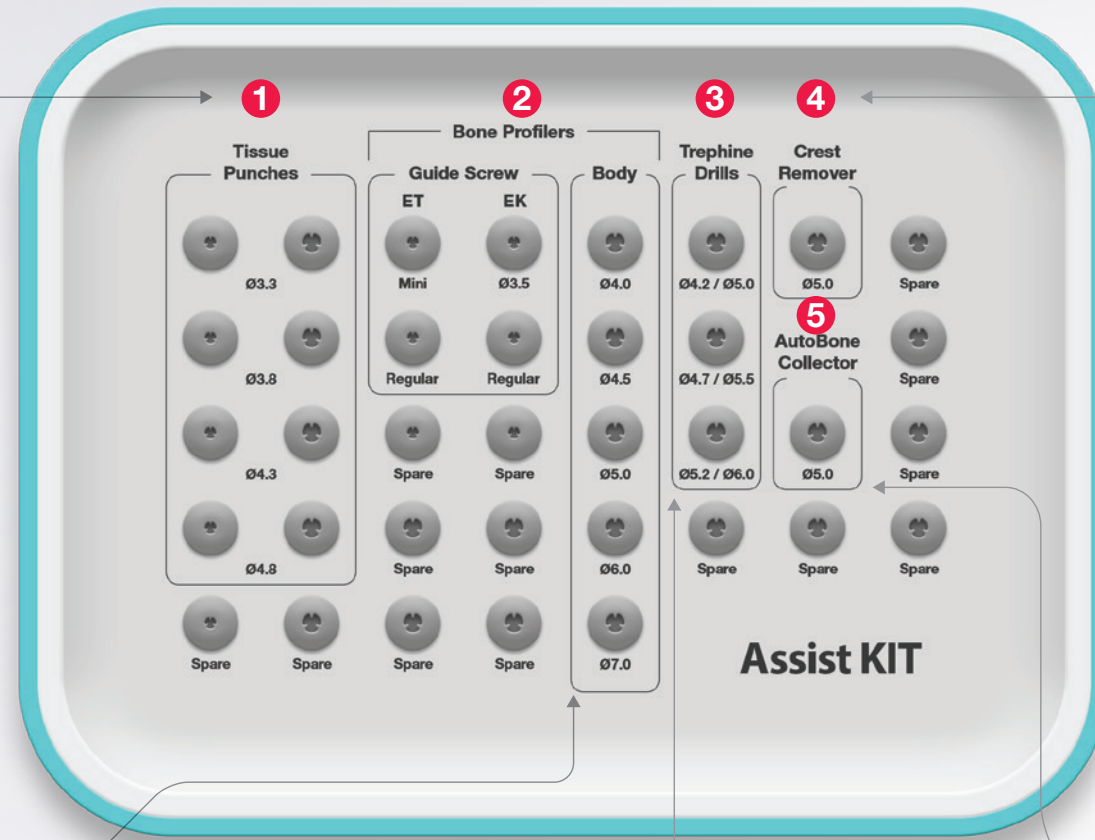

HIOSSSEN

IMPLANT

Smiles that last a lifetime






Please contact your local sales representative or visit our website today to learn more about Hiossen and its products.



All Hiossen Implants are processed and
Manufactured in the USA

270 Sylvan Ave. Ste 1130, Englewood Cliffs, NJ 07632
Marketing@hiossen.com | 888.678.0001 | www.hiossen.com

 @Hiossen  @HiossenImplants  @Hiossen

Copyright © 2024 Hiossen Inc. All rights reserved.

All information included herein is subject to change without notice. No part of this brochure may be reproduced in any manner without the prior written permission of Hiossen Inc. Hiossen and the Hiossen logo are trademarks of Hiossen Inc. All other brand or product names are trademarks of their respective companies or organizations.

PF24ASKLTR1.1