



Meet Hiossen
ET Implant System

HIOSSEN
IMPLANT

ET (Essential Treatment) Implant System

Hiossen ET stands for **Essential Treatment**, as it provides all the necessary tools that a clinician needs in order to perform a successful treatment.

Our ET System has become one of the best choices for dental implant systems in the industry due to our commitment to extensive research and the development of high quality, safe, and innovative implants.

Immediate results

Hiossen ET Implants design has a tapered body providing outstanding primary stability and surgical simplicity.

Hiossen ET System



Excellent primary stability

ET Implant system has been specifically designed to achieve predictable and minimally invasive placements in all bone types



Clinical data and ongoing clinical trials

Applying strict quality standards based on technology and regulations that are needed by country and by product



Prosthetic portfolio

The scope of our prosthetics covers a broad range of patients' conditions. By using our original components, you can have peace of mind to expect perfect matching and high-end manufacturing precision. Our SmartFit abutments helps you achieve a high degree of aesthetic needs and create an ideal smile design



Hiossen surface treatment

- **SA Surface:** Sand-blasted with Alumina and Acid-etched surface treatment

Hiossen ET features


Our Implants stand for maximized clinical performance, improved esthetics and high patient satisfaction.

Our world-class implant surface treatment SA helps by providing a minimally invasive procedure, excellent predictability, and treatment success.

ET III

1.5 Degree taper body

The apex mimics natural tooth providing strong initial stability even for immediate cases

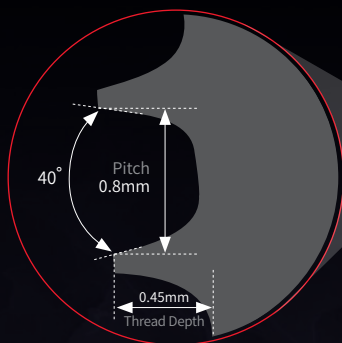

Optimal for all bone types

Suitable for immediate and early loading


Primary usage

Any type of cases.
Ultra wide implants can also be used for rescue cases

ET IV

6 Degree taper body

Minimized the damage to adjacent teeth roots. Similar shape to natural tooth root

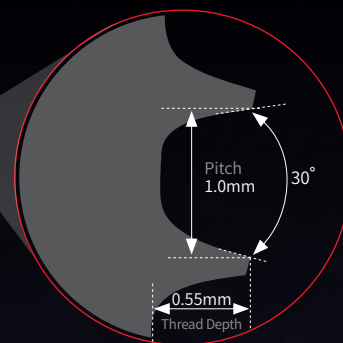

*Based on ETIV Ø4.5

Optimal for soft bone

In order to insert the implant into the osteotomy, higher torque is required when the implant body design is more tapered.

ETIV's taper body has a wider pitch than ETIII's to minimize over-torque.

Additional specifications



Pitch for ETIV Ø4.0 : 0.8mm
Pitch for ETIV Ø4.5, Ø6.0, Ø7.0: 1.0mm
Pitch for ETIV Ø5.0: 1.2 mm

* Thread depth varies with each size

ET Implant system macro design contributes to the process of Osseointegration

Submerged internal taper connection


Internal hex submerged Implant.
11° morse taper design facilitates conical sealing and stabilizes the implant more so that the chances of screw loosening is minimal.
One-stage and two-stage surgery are applicable.

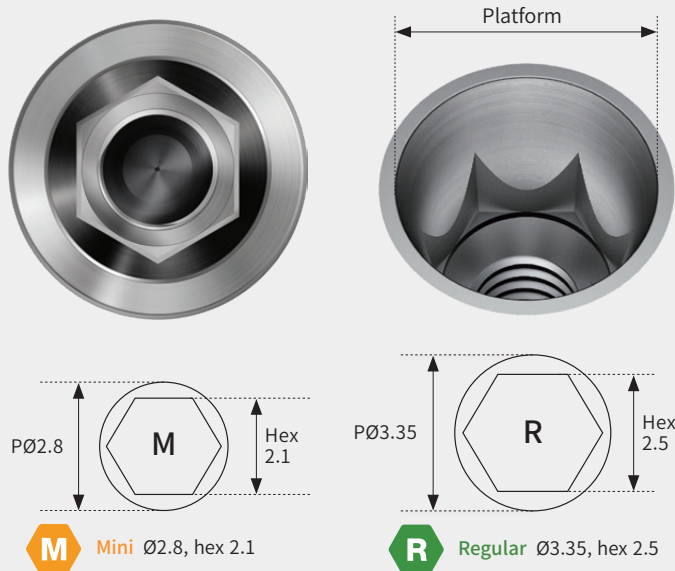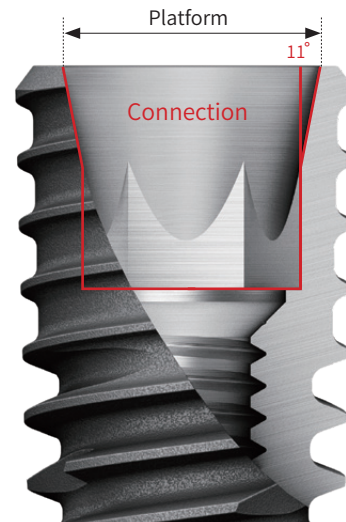





Strong internal hex connection

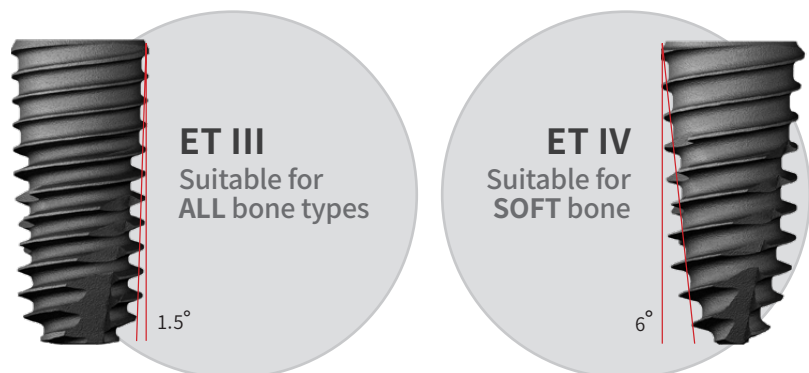The Hex connection distributes mastication force to both abutment and implant evenly and prevent micro-movement while increasing stability.




Two options are available depending on the clinical indications



Specifications for **ET** Implant System




CONNECTION	PLACEMENT	DELIVERY OPTIONS	BODY DIAMETERS
 <p>Internal hex conical connection</p>	 <p>1.0mm Sub-crestal</p>	 <p>No Mount Pre-Mounted</p>	 <p>Diameter</p> <ul style="list-style-type: none"> M Mini: 3.2mm, 3.5mm R Regular: 4.0mm, 4.5mm, 5.0mm, 5.5mm R Ultra Wide: 6.0mm, 7.0mm

Specifications for **ETIII**

BODY TYPE	HEX PLATFORM	APICAL DIAMETER	LENGTHS (mm)
 <p>1.5° Tapered Body</p>	<ul style="list-style-type: none"> M Mini: 2.1mm (Narrow) R Regular: 2.5mm 	 <p>Apical Diameter</p> <ul style="list-style-type: none"> M Mini: 2.2mm, 2.5mm R Regular: 2.8mm, 3.1mm, 3.7mm, 4.0mm, 4.2mm, 5.1mm 	 <p>Length</p> <ul style="list-style-type: none"> M Mini: 8.5 10 11.5 13 15 R Regular: 7 8.5 10 11.5 13 15 R Ultra Wide: 6 7 8.5 10 11.5 13

**Note: For Implants with a length of 6mm are only available for implant with Ø5.0 above.
For Implants with a length of 15mm can only be placed as custom orders.

Specifications for **ETIV**

BODY TYPE	HEX PLATFORM	APICAL DIAMETER	LENGTHS (mm)
 <p>6° Tapered Body</p>	<ul style="list-style-type: none"> R Regular: 2.5mm 	 <p>Apical Diameter</p> <ul style="list-style-type: none"> R Regular: 1.8mm, 2.0mm, 2.2mm, 2.9mm, 3.8mm 	 <p>Length</p> <ul style="list-style-type: none"> R Regular: 7 8.5 10 11.5 13 R Ultra Wide: 7 8.5 10 11.5 13

**Note: Product availability varies by country depending on the approval status by the regulatory authority for medical products in each country

Provide a complete tooth replacement to your patients

Our extensive portfolio conveniently offers a complete range of clinical indications and choices of your need.

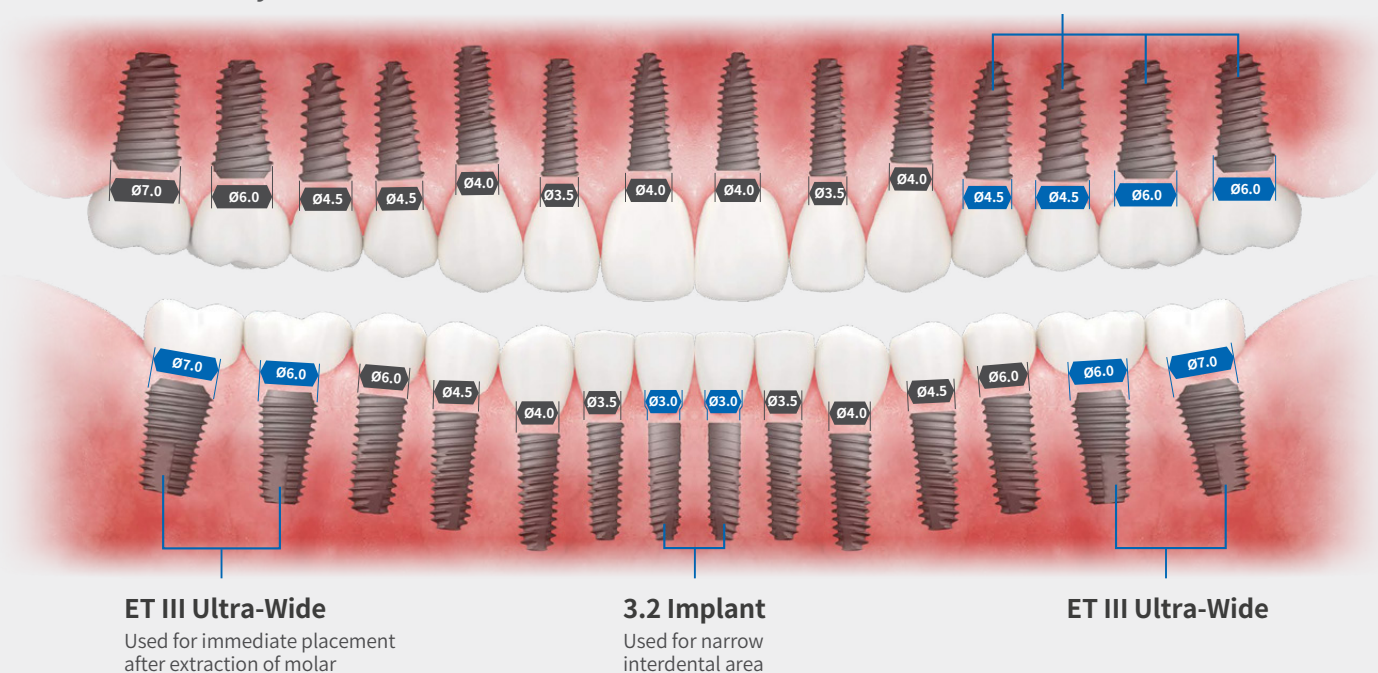


Fig 1. Implant size selection. (Hiossen's recommendations)

Actual clinical conditions and clinician's assessment of the patient should be the main criteria for choosing the size of an implant for a particular area.



Hiossen Mini Narrow 3.2 Implant

The ET 3.2 is made of Ti-alloy

- Specifically designed for restrictive spaces of the esthetic zone, like the lateral incisors.
- Offers a wider restorative options and ideal esthetic results.
- Can be used to mount full or partial dentures, crowns and other restorations.



Hiossen Mini 3.5 Implant

The ET 3.5 Mini Implant is made of "Grade 4" Unalloyed Titanium

- Less impact on the jawbone and requires less bone for support.
- Can be used to mount full or partial dentures, crowns and other restorations.

****Note:** Product availability varies by country depending on the approval status by the regulatory authority for medical products in each country

Optimized predictability and treatment success

Essential treatment system line-up.

Careful selection of the implant is mandatory to achieve successful, predictable, and a consistent treatment outcomes.

IMPLANT TYPE	IMAGE	IMPLANT FEATURES
ET III		<p>1.5° tapered body angle is ideal for obtaining primary stability.</p> <p>Suitable for immediate and early loading.</p> <p>Suitable for all types of bones.</p>
ET IV		<p>ETIV is designed with more tapered body and with wider pitch to achieve stronger primary stability while minimizing the loss of implant stability at compression site.</p> <p>Suitable for soft bone and maxillary sinus procedures.</p>



Hiossen Regular Implant

The ET Regular Implants are made of “Grade 4” Unalloyed Titanium

- Provides more primary stability.
- Permanent solution to tooth loss.
- Recommended for all bone types and procedures.



Hiossen Ultra-wide Implant

The ET Ultra-Wide Implants are made of “Grade 4” Unalloyed Titanium

- Alternative when the bone quality in the posterior teeth is relatively low.
- Rescue cases.
- Suitable for posterior, mainly molars.

****Note:** Product availability varies by country depending on the approval status by the regulatory authority for medical products in each country

Hiossen IMPLANT

Smiles that last a lifetime

Please contact your local sales representative or
visit our website today to learn more about Hiossen and its products.



All Hiossen Implants are processed and
Manufactured in the USA

270 Sylvan Ave. Ste 1130, Englewood Cliffs, NJ 07632
Marketing@hiossen.com | 888.678.0001 | www.hiossen.com

@hiossen_implants

@HiossenImplants

@Hiossen

Copyright © 2022 Hiossen Inc. All rights reserved.

All information included herein is subject to change without notice. No part of this brochure may be reproduced in any manner without the prior written permission of Hiossen Inc. Hiossen and the Hiossen logo are trademarks of Hiossen Inc. All other brand or product names are trademarks of their respective companies or organizations.

PB22ETSLTR1.0