



Easy Maintenance

- Easy disinfection with simplified motor separation

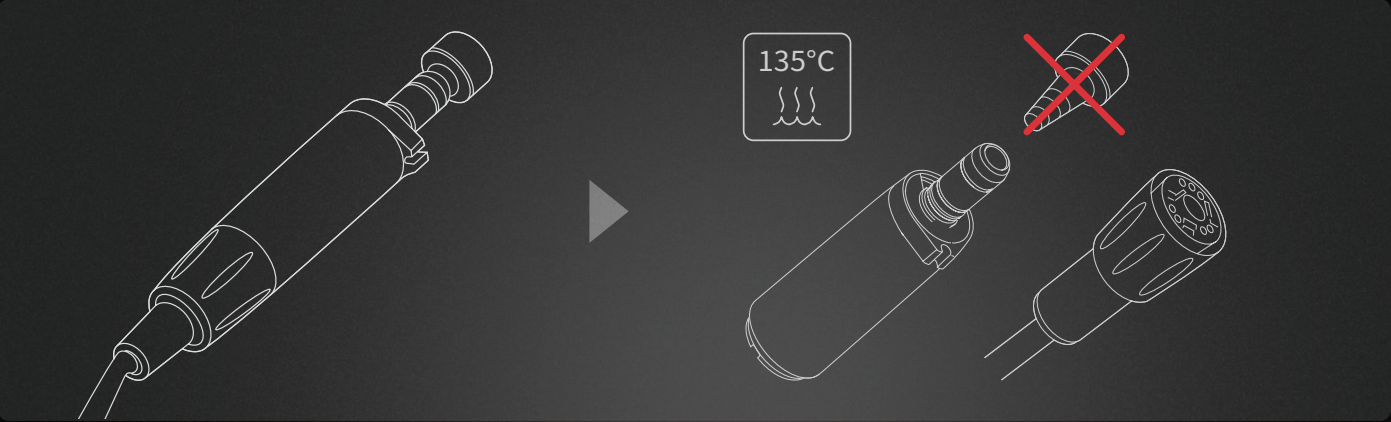

- Easy cleaning with an optimized removable foot pedal

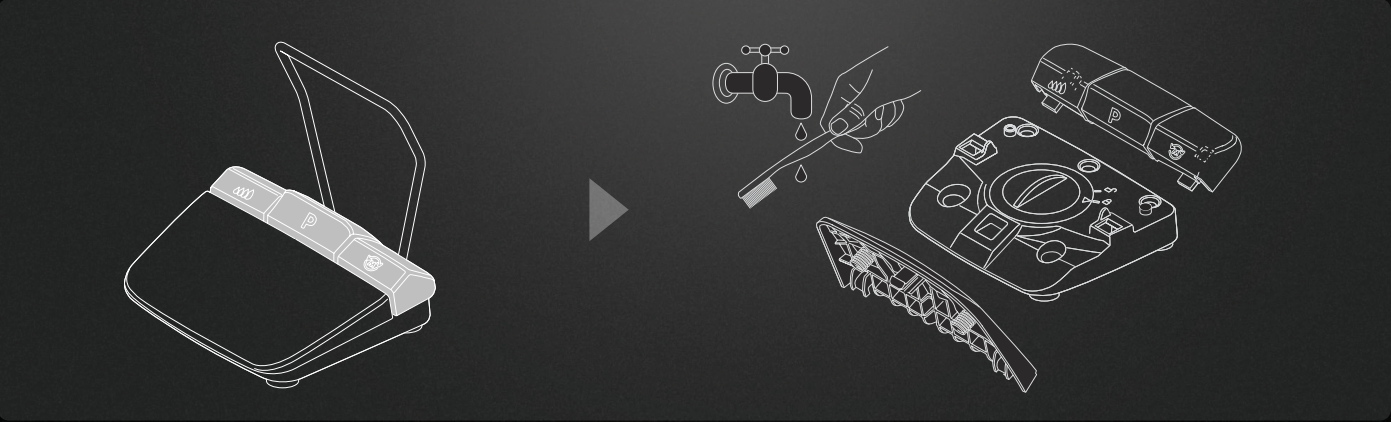

The best performance can be maintained when used with the genuine KaVo Oil Spray

Technical Specification

Length x Width x Height	265 x 255 x 100mm
Weight	1.9kg (Foot pedal 1.1kg / Motor 125g)
Input voltage	100 ~ 240v
RPM speed	300 ~ 40,000rpm
Motor max. torque	80Ncm
Pump speed	30 ~ 110ml/min
Cable length	Foot pedal 2.5m / Motor 2m

The Best Choice for Optimal Care

SM5

- Enhanced power output
- More compact and lightweight
- Constant torque regardless of the RPM




Enhanced Engine Power

- More powerful compared to INTRAsurg 300 Plus
- Smaller and lighter engine, minimizing wrist fatigue during surgery




Intelligent Motor Control

- The SM5 motor provides consistent torque no matter the RPM, ensuring greater safety and consistency during implant insertion

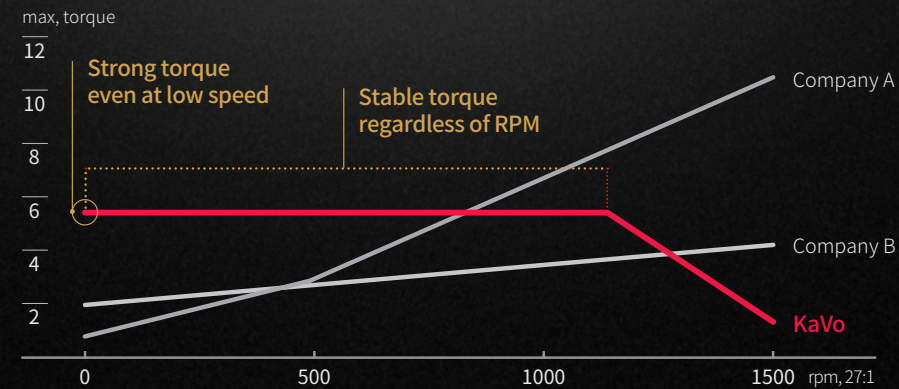


Flexible Tubing

- Thinner tubing for a more flexible operation
- Enhanced durability of the tubing joint

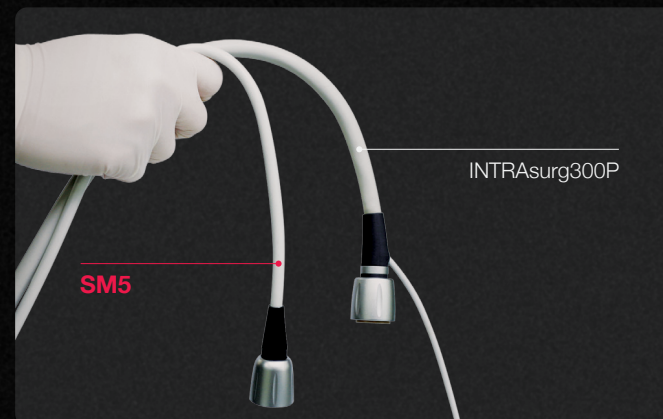

One-touch Auto Calibration

- When changing the engine torque value, the new setting can be calibrated with a single touch of a button without any additional equipment

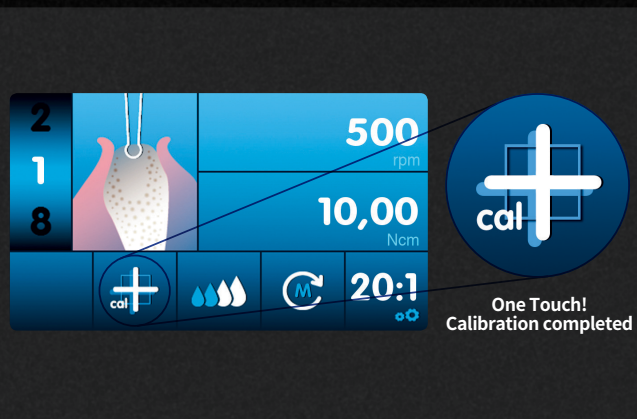

Surgical Instruments

Contra-angle 20:1



S201L* LED
S201C

S201XL* LED
S201XC

Straight-angle 1:1



S11L LED
S11C

- Optimized clamping system - avoid interfering with tool shafts
- EXPERTgrip - improved grip, comfort and hygiene
- Easy maintenance with a removable head
- External and internal irrigation available

*Optical fiber: reliable viewing of the surgical field

Max. speed : 15 ~ 2,000rpm
※ Basic Components : S201L

- Can be used with a standard length 44.5mm handpiece bur and contra-angle bur
- Used with external irrigation

Max. speed : 40,000rpm