

Discover the innovation that sets it apart.

Let us handle the rest, while you prioritize what truly counts.



**Free** Access to  
**Medit Link & Medit Apps**




**3-Year** Warranty  
(Applies Only to i900)



User-Friendly  
**Medit Community**  
(Facebook User-Page & Medit Academy)

Specification

Category		Medit i900
Scanning Technology	Imaging Technology	3D-in-motion video technology   3D full color streaming capture
Features	Field of View (FOV)	Large: 18 x 15 mm   Mid: 14 x 11.5 mm   Small: 10 x 8.2 mm
	Touch Interface	Yes
	Remote Control Mode	Yes
	UV-C Disinfection	Yes
	Unibody	Yes
	Anti-Fogging Technology	Air Flow
Handpiece	Dimensions (W x H x L)	36.7 x 35.3 x 223.4 mm
	Total Weight (Incl. Tip)	165 g
<div>Tip</div>	Size	Full Tip (W x H x L)
		Tip Head (W x H)
Autoclavability	Autoclave Cycles	Up to 150 times
	Autoclave Type	Sterilize for 4 mins at 134°C and dry for 20 min (Pre-Vacuum Type) Sterilize for 10 mins at 135°C and dry for 30 min (Gravity Type)
Connectivity		Wired
Cable		USB C-Type Power Delivery Cable: 2m & 2.5m

For more information or to request a demo visit us at <https://www.medit.com/contact-us/>

# MEDIT i900

## The Ultimate Solution, That Does It All!

Our premium scanner **Medit i900** redefines the scanning process, making it not just a task but a **fascinating journey from purchase to performance** for every user.





## 360° Touch Band

- Double-Tap to Start & Stop Scan
- Swipe to Navigate Between Scanning Phases



## 3<sup>rd</sup> Gen New Optical Engine

- Improved Metal Scanning
- Enhanced Scan Accuracy
- Deeper Scanning Depth
- Better Texture



*Delivers Accurate Scan Bodies with Ease*



*Provides Precise Inlay Preparation*

## Unibody

Seamless Handpiece



ALL IN for Your New Scanning Experience

# MEDIT i900



## Hygienic Touch Pad

No Physical Buttons or Gaps to Reduce the Risk of Cross-Contamination



## Our Lightest, Most Compact

Designed for Everyone  
on Your Dental Team



## Vibrant LED Indicators

Stay Informed From Connection to Scanning

Ready for Scan

Scanning

UV-C LED in Process


Caution

Alert




## Wider Field of View 18 x 15mm (552 x 460px)

Get Precise Scans Faster with Fewer Image Captures

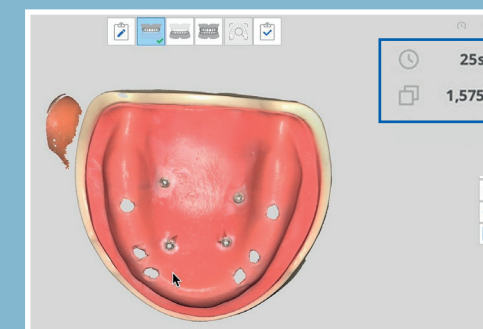

*Medit i900*

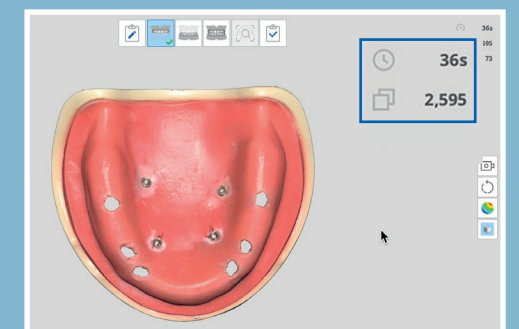

*Comparative Product*