



Performance You Can Trust, Innovation You Can Feel


K5 Unit Chair

- The K5 unit chair features cutting-edge technology, including intuitive controls and advanced water sterilization system
- Engineered to withstand the demands of daily use, ensuring longevity and reliability
- Sleek and modern design enhances the overall look of your dental practice

HIOSSSEN
IMPLANT


Slim Backrest

- Ensures an ideal treatment position by reducing obstruction to the dentist's knee

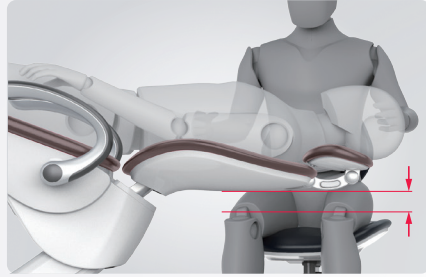

Seat Tilting Function

- Maintains a stable patient treatment position when the backrest reclines

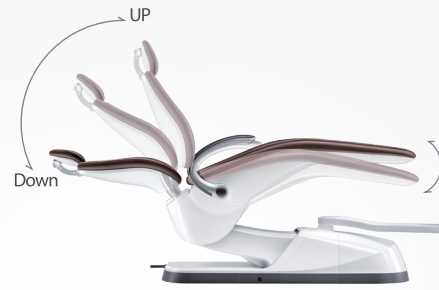

Comfortable Seat

- Provides excellent support for both the patient and the dentist, helping to reduce strain during procedures



Two-Jointed Head Rest



- The angled headrest is designed to support children, the elderly, and patients in wheelchairs

Doctor's Stool



- The backrest can be tilted to provide adjustable support for the dentist's hips during treatment
- The seat base can be tilted forward to reduce pressure on the thighs, minimizing leg fatigue

HANARO Console

- The integrated design offers convenient storage for personal items, paper cups, and tissues, maximizing space efficiency


Rotatable Cuspidor

- The cuspidor can be adjusted to allow the patient to comfortably rinse their mouth after treatment



Joystick

- Multi chair position adjustment

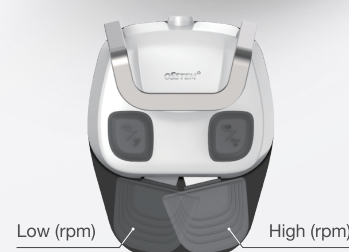


LED Light & Sterilization

- Blue light Cut-Off Function: Controls the resin curing time
- On/Off sensor for hands-free light control
- Autoclavable Handles



Adjustable Hand-piece Cradle

- Protects the dentist's hand by reducing interference with the surgical instruments

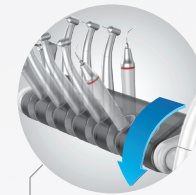

Warm Water Supply System

- Provides temperature control on the hand-piece spray



Doctor's Table

- Left handed option available

Wireless Foot Controller

- Handpiece RPM can be easily adjusted by turning the foot pedal left or right as needed


Autoclavable Doctor's Table-Handle

- Can be detached and sterilized for easy cleaning



Sterilized Water Supply Device

- Supplies sterilized water to effectively suppress bacteria and biofilm in the water pipe

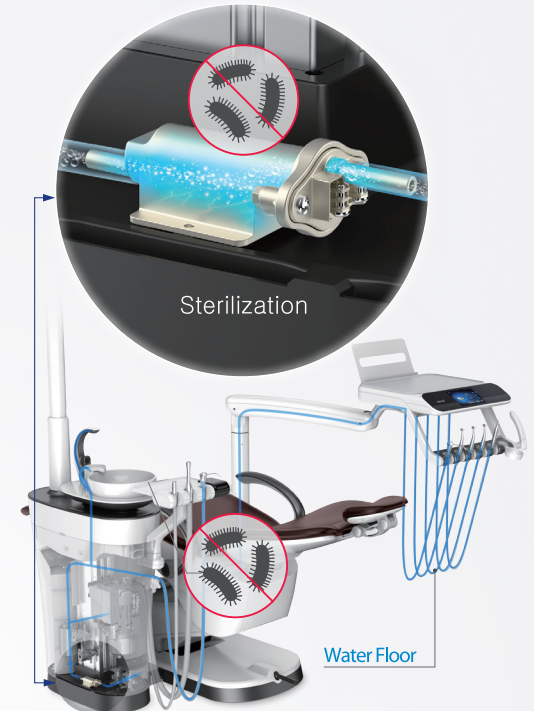

Built-in Water Pipe Cleaning System




- Oxygenated water assists in effectively cleaning the water pipe



K5 Unit Chair Instruments

	Mini Head	Standard Head		Coupling
High-speed Handpieces	NSK H511L - M900L	H311L - M800L		PTL-CL
	Bien Air Tornado SLK	Tornado LK	Bora 2 LED	Bora 2 LK

	Contra Angle Handpiece		Straight Handpiece
Low-speed Handpieces	NSK E200IL	E205L 1:5	E400i
	Bien Air CA 1:1L	CA 1:5L	HP 1:1L
Electric Motor	EIL (Bien-Air)		

Scalers			
			
	Dmetec Scaler	Osstem S5 (EMS)	Osstem S5L (EMS)

Product Details

Color Options



Specifications

Power Supply Voltage	AC 110~120 V, 60 Hz
Power Consumption	1.2 KVA
Seat Elevation Method	Vibrational
Seat Height	Max. 2.8 ft, Min. 1.5 ft
Backrest Angle	Maximum angle 71°
Weight	441 lbs
Maximum Load	331 Lbs
Monitor Size	Smart 27"

Actual Dimensions

Width	5.0 ft
Length	6.1 ft
Height	6.3 ft

