

K5 Unit Chair Instruments


Air Handpieces	Standard Head	Mini Head		Coupling
	NSK SMAX 900L	SMAX 800L		PTL-CL LED
	Bien Air Bora 2 LK			UNIFLEX

Electric Handpieces	Contra Angle Handpiece		Straight Handpiece
	NSK SMAX M25L	SMAX M95L	SMAX M65
	Bien Air Classic CA 1:1L	Classic CA 1:5L	Classic PM1L
	Electric Motor Bien-Air MCX		

Scalers	Dmetec Scaler
---------	---------------

Product Details

Recommened Room Dimensions: 8 ft x 10 ft

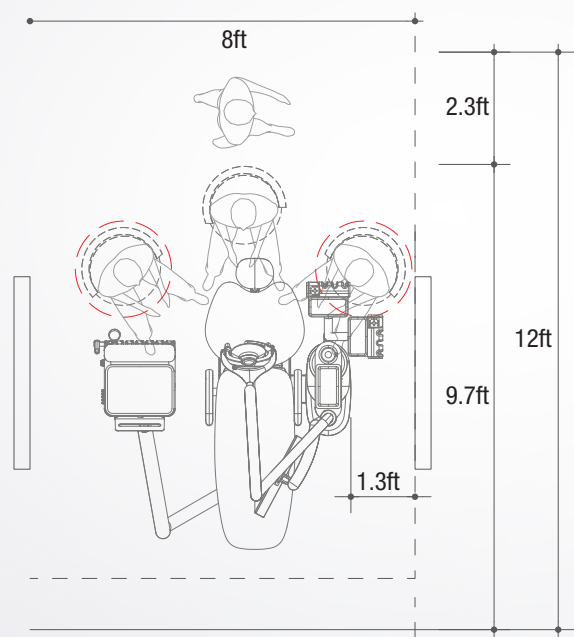

Color Options



Specifications

Power Supply Voltage	AC 110~120 V, 60 Hz
Power Consumption	1.2 KVA
Seat Elevation Method	Vibrational
Seat Height	Max. 2.8 ft, Min. 1.5 ft
Backrest Angle	Maximum angle 71°
Weight	441 lbs
Maximum Load	331 Lbs
Monitor Size	Smart 27"




Performance You Can Trust, Innovation You Can Feel

K5 Unit Chair

- The K5 unit chair features cutting-edge technology, including intuitive controls and advanced water sterilization system
- Engineered to withstand the demands of daily use, ensuring longevity and reliability
- Sleek and modern design enhances the overall look of your dental practice

Slim Backrest

- Ensures an ideal treatment position by reducing obstruction to the dentist's knee

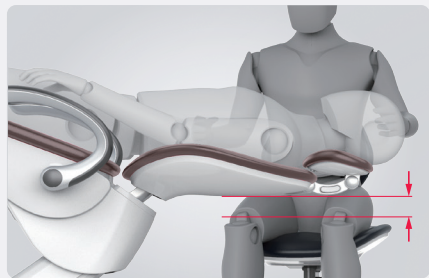

Comfortable Seat

- Provides excellent support for both the patient and the dentist, helping to reduce strain during procedures



Two-Jointed Head Rest



- The angled headrest is designed to support children, the elderly, and patients in wheelchairs


Doctor's Stool



- The backrest can be tilted to provide adjustable support for the dentist's hips during treatment
- The seat base can be tilted forward to reduce pressure on the thighs, minimizing leg fatigue

Seat Tilting Function

- Maintains a stable patient treatment position when the backrest reclines

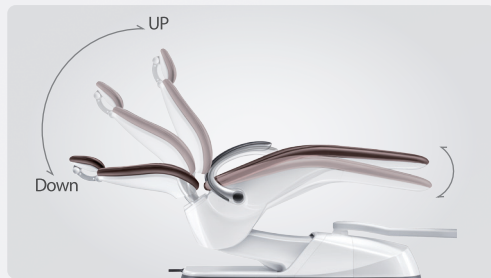

Rotatable Cuspidor

- The cuspidor can be adjusted to allow the patient to comfortably rinse their mouth after treatment



HANARO Console

- The integrated design offers convenient storage for personal items, paper cups, and tissues, maximizing space efficiency



LED Light & Sterilization

- Blue light Cut-Off Function: Controls the resin curing time
- On/Off sensor for hands-free light control
- Autoclavable Handles




Autoclavable Doctor's Table-Handle

- Can be detached and sterilized for easy cleaning



Adjustable Hand-piece Cradle

- Protects the dentist's hand by reducing interference with the surgical instruments

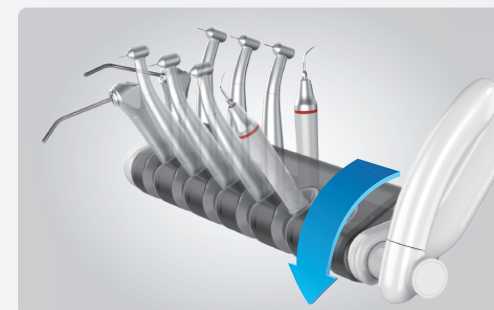

Doctor's Table

- Left handed option available



Warm Water Supply System

- Provides temperature control



Built-in Water Pipe Cleaning System (Optional)

- Oxygenated water assists in effectively cleaning the water pipe



Joystick

- Multi chair position adjustment



Low (rpm) High (rpm)

Wireless Foot Controller

- Handpiece RPM can be easily adjusted by turning the foot pedal left or right as needed