



The Best Choice for Optimal Care


# SM5

- Enhanced power output
- More compact and lightweight
- Constant torque regardless of the RPM

**HIOSSEN**  
IMPLANT


## Enhanced Engine Power

- More powerful compared to INTRASurg 300 Plus
- Smaller and lighter engine, minimizing wrist fatigue during surgery

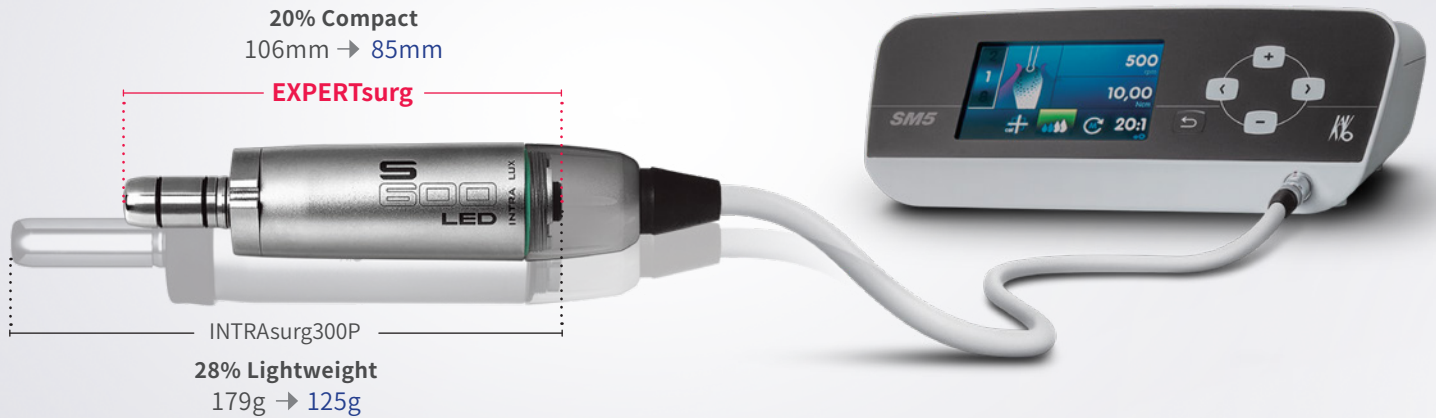

## Intelligent Motor Control

- The SM5 motor provides consistent torque no matter the RPM, ensuring greater safety and consistency during implant insertion

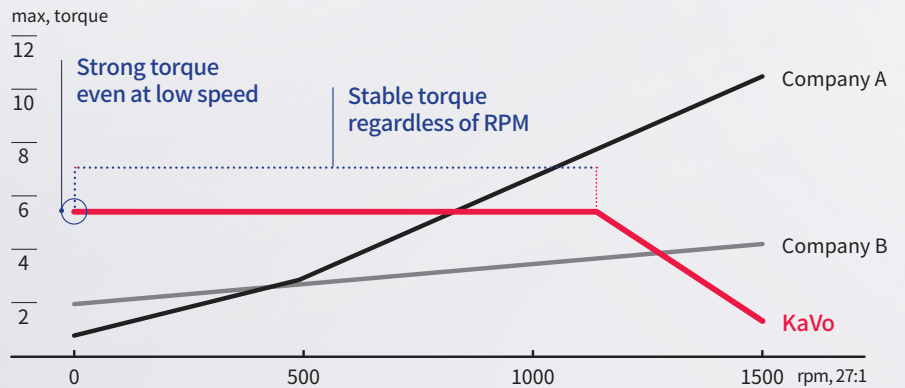


## Flexible Tubing

- Thinner tubing for a more flexible operation
- Enhanced durability of the tubing joint



## One-touch Auto Calibration

- When changing the engine torque value, the new setting can be calibrated with a single touch of a button without any additional equipment

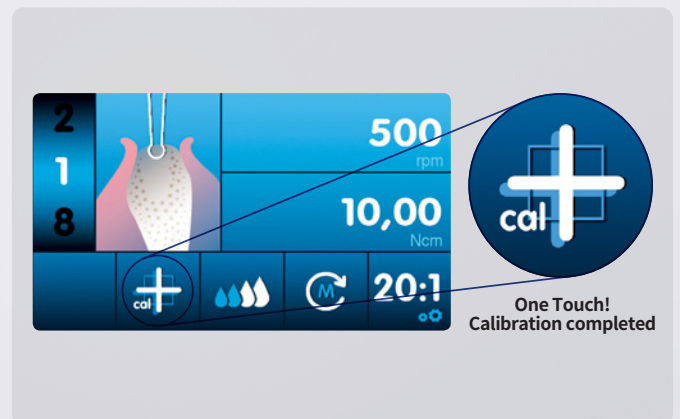

## Surgical Instruments

Contra-angle 20:1

---



S201L\* LED  
S201C

S201XL\* LED  
S201XC

- 
- Optimized clamping system - avoid interfering with tool shafts
  - EXPERTgrip - improved grip, comfort and hygiene
  - Easy maintenance with a removable head
  - External and internal irrigation available

\*Optical fiber: reliable viewing of the surgical field

**Max. speed : 15 ~ 2,000rpm**

※ Basic Components : S201L

Straight-angle 1:1

---



S11L LED  
S11C


- 
- Can be used with a standard length 44.5mm handpiece bur and contra-angle bur
  - Used with external irrigation

**Max. speed : 40,000rpm**



## Easy Maintenance

- Easy disinfection with simplified motor separation

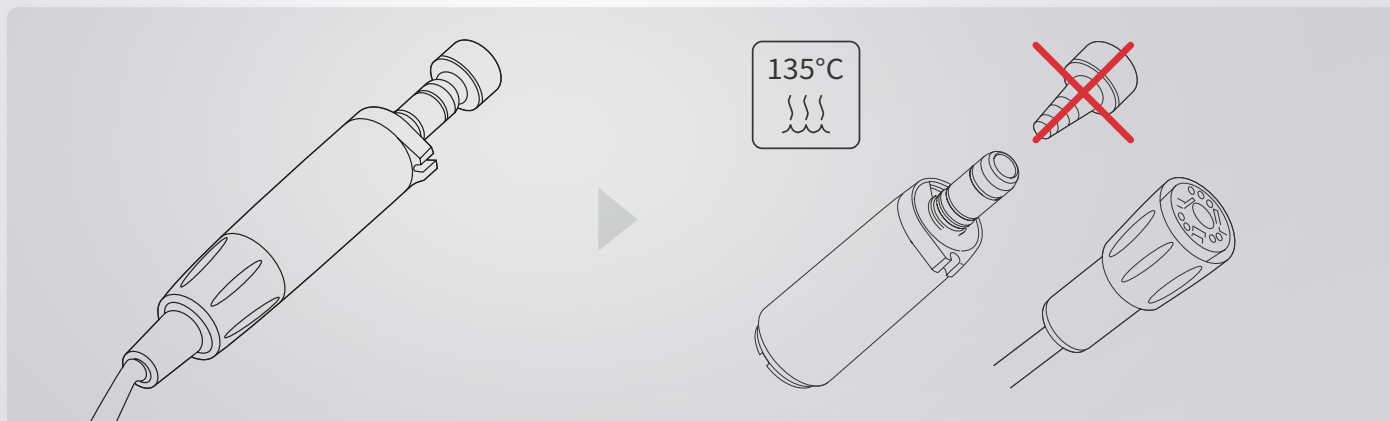

- Easy cleaning with an optimized removable foot pedal



The best performance can be maintained when used with the genuine KaVo Oil Spray

## Technical Specification

Length x Width x Height	265 x 255 x 100mm
Weight	1.9kg (Foot pedal 1.1kg / Motor 125g)
Input voltage	100 ~ 240v
RPM speed	300 ~ 40,000rpm
Motor max. torque	80Ncm
Pump speed	30 ~ 110ml/min
Cable length	Foot pedal 2.5m / Motor 2m