

Five Essential Tools  
to Achieve Surgical Convenience

# Assist Kit

- 1. Tissue Punch
- 2. Bone Profiler
- 3. Trephine Drill
- 4. Crest Remover
- 5. Autobone Collector



**HIOSSEN**  
IMPLANT

## 1 Tissue Punch



### Minimizes incisions and reduces healing time

- The guide pin configuration enables precise tissue removal at the designated location
- Recommend using a diameter 0.7-1.5mm smaller than the healing abutment diameter




## 2 Bone Profiler



### Removes protruding bone without internal damage of the implant

- Easily removes adjacent obstructive bone around the implant
- Detachable guide protects interior of the implant while removing the bones


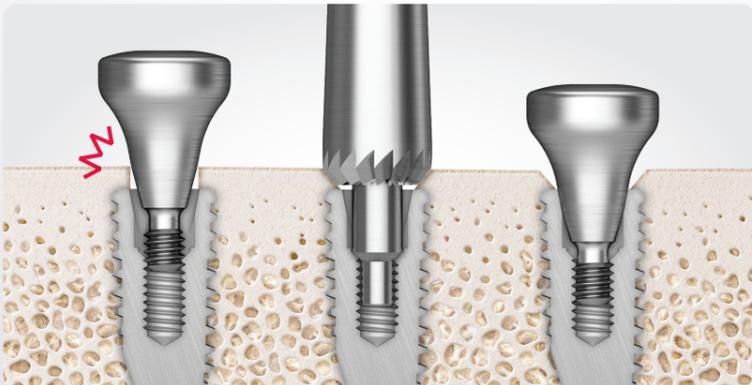


## 3 Trepine Drill



### Thin blade wall (0.4mm) minimizes amount of the bone removed

- The drill is designed to harvest larger block bone from the donor site or to remove the failed implant
- The 0.4mm thin wall design reduces the quantity of removed bone

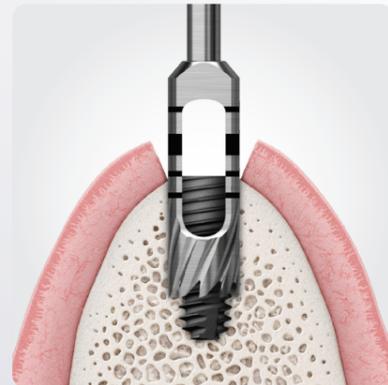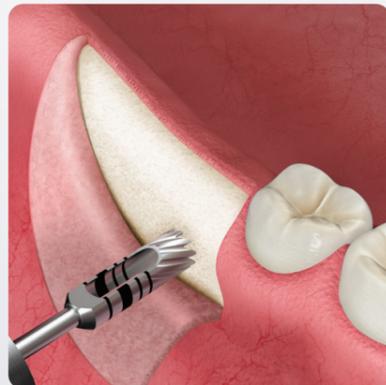


## 4 Crest Remover



### Net-shaped blades of the bur speeds up ridge widening

- The blade is designed to minimize bone splattering effectively during the drilling process

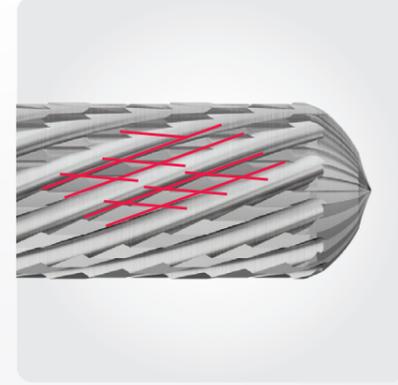


## 5 Autobone Collector



### Collects bone chips through drilling

- Efficiently gather autogenous bone chips with a strong cutting force and design
- The stopper enables a safer method for collecting bone material (0.18g acquired at one operation)

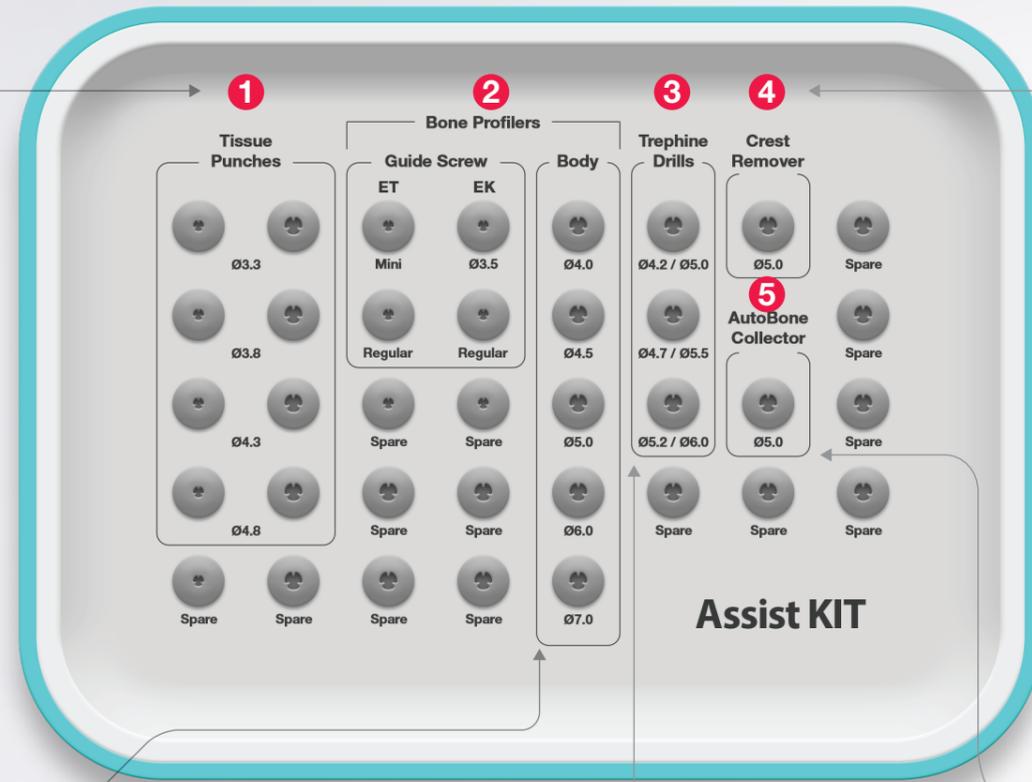

# **Hiossen** IMPLANT

Smiles that last a lifetime



Please contact your local sales representative or visit our website today to learn more about Hiossen and its products.



All Hiossen Implants are processed and  
Manufactured in the USA

105 Challenger Rd. 7th Floor, Ridgefield Park, NJ 07660  
Marketing@hiossen.com | 888.678.0001 | www.hiossen.com

 @Hiossen       @HiossenImplants       @Hiossen

Copyright © 2026 Hiossen Inc. All rights reserved.

All information included herein is subject to change without notice. No part of this brochure may be reproduced in any manner without the prior written permission of Hiossen Inc. Hiossen and the Hiossen logo are trademarks of Hiossen Inc. All other brand or product names are trademarks of their respective companies or organizations.

PF26ASKLTR1.0